
ZebraScope® 5 Limited-Reusable Ureteroscope



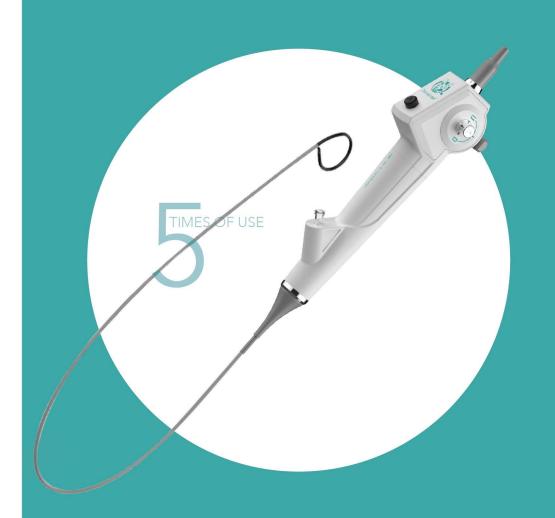

THE LIMTED-REUSABLE DIGITAL URETOROSCOPE 5 TIMES USE

KEY BENEFITS

- Up to 5 times sterilizations
- Maintaining the best performance and imaging
- Reducing the cost of disposable ureteroscope adoption
- Usage timer available

ZebraScope® 5 limited-reusable Uretoroscope is a disposable uretoroscope with 5 times limited use solution that gives you a way to take control of your clincal budget. This is to meet the needs of local hospitals which cannot afford the traditional expensive reusable ureteroscope and disposable ureteroscope.

As the first brand of single-use flexible ureteroscope launched in China, the successful launch of ZebraScope® is available in over 2000 hospitals. With the awareness of coming from the clinic and serving the clinic, we discovered that there are many hospitals and clinics that lack patients due to a lack of equipment worldwide. In Brazil, according to RH Astolfi (06, 2020), the main reason discussed that the local hospitals were hugely affected due to a lack of equipment or adequate hospital infrastructure with for URS (23%) and F-URS (38.1%), while for PCNL the lack of practical experience was the main reason (29.9%). Carrying the mission of being the provider of the disposable endoscope to improve human health, we launch the ZebraScope® 5 limited reusable ureteroscope which combines the features of both disposable and reusable ureteroscope.

ZebraScope® 5 limited-reusable ureteroscope enables hospitals to offer kidney stone management locally with an affordable budget. It can be used up to 5 times during sterilizations while maintaining the

best performance conditon of the scope. As a result, ZebraScope® 5 simplified workflow, provides high-quality images and reduces cross-contamination, so you can assure your patients get a local kidney treatment with the best performance.

Environmental

Happiness Works Medical is committed to a healthy environment. With this in mind, we have designed the ZebraScope® and the packaging to reduce the environmental impact. We have also kept the environment in mind when designing our packaging which is minimized to create less waste. The cardboard and plastic packaging is 100% recyclable.

SPECIFICATIONS

OPTICAL SYSTEM		
Field of View	100°	
Direction of View	0°	
Depth of Filed	2mm - 50mm	
Illumination Method	thod LED	
INSERTION PORTION		
Bending Angle	275° up l 275° down	
Insertion Cord Diameter	≤2.9mm	
Distal End Diameter	7.4 Fr	
Length of Insertion	670mm±10mm	
WORKING CHANNEL		
Minimum Instrument Channel Diameter	3.6 Fr	
STERILIZATION		
Method of Sterilization	EO	
MATERIAL		
DEPH Phthalates	Does not contain any DEPH phtalates	
No allergenic substances	Not made with natural rubber latex	
DIMENSIONS		
Weight	298g	

VIDEO PROCESSOR HV-200					
Voltage	AC 100-240V				
Frequency	50/60Hz	Processor Front			
HQ Pixel	200*200,400*400,800*800,1000*1000	Happiness vanks **A.Si medical** W-200			
Key Function	LED Dimming Adjustment	Host Switch			
	Image zooming & freezing				
	Dislay Adjustment Brightness & Saturation & Contrast	Remo Control panel Socket			
	White Balance Adjustment	Socket			
	Picture & Video Taking				
	Video Ratio	Processor Back			
	Handle Bunton Customization				
Video Output	Support digital signal DVI, HDMI, SDI output	Grounding			
	Pressure measurement mode output	POWER -200 SOFE USB RS202 DVI HIMI SV 💠 POST			
Dimensions	300mm X 300mm X 100mm	Power Video Output Interface interface			
Weight	3.109kg				

ORDER SPECIFICATIONS

PRODUCT	TYPE NO.	QUANTITY
ZebraScope® 5 Limited-Reusable Ureteroscope	HUR-101	1pcs
ZebraScope® Processor	HV-200	1pcs



Shortcut Button

The shortcut button is located on the top of the handle for easy operation. During the procedure, you can quickly shoot, video, freeze the frame, zoom out, and white balance



The only one equipped with snake-bone-covered woven mesh

ZebraScope is the only endoscope containing snake bone and covered with woven mesh which is safer for clinical use. It would effectively prevent accidental damage caused by the broken snake bone which could lead to the risk of being unable to withdraw the ureteroscope.



Unique stereo deflection function

The distal end with a unique stereo deflection function of 12.5° and a minimum bending radius of 8mm which enables the ureteroscope with the best performance on the small IPA renal calyx stones.



Usage TimerUsage timer available



Bending Radius 8mm



Passive bending



Deflection 275° up 275° down



Stereo deflection angle 12.5°



Unique snake bone design Covered with woven mesh



Anhui Happiness Workshop Medical Instruments Co., Ltd

Building 2, Wisdom Industrial Park, Bengbu Jiangqiao Road, Bengbu, Anhui, China 233000 Specifications indicated in the above table are subject to change without any notice.
All trademarks are the property of their respc tive owners.

www.hwmeds.com hello@hweds.com